

healthLINK^{SIX}6

SECURE INFORMATION EXCHANGE

BETTER PRACTICE - BETTER CARE

Installation Guide Companion

Applicable to setup Impulse for PIT messages

Healthlink Messaging Systems (HMS) 6.6.x



© HealthLink 2011. All rights reserved. No reproduction, transmission, transcription, storage in a retrieval system, or translation into any language or by any means, electronic, mechanical, optical, chemical, manual, or otherwise, any part of this document without express written permission of HealthLink Limited.

HealthLink, the HealthLink logo, HealthLinkSIX, HealthLinkSIX Messaging System (HMS) are trademarks of HealthLink Limited in New Zealand and Australia. All other brand names are trademarks or registered trademarks of their respective companies

Liability Notice: Every effort has been made to ensure that the information in this document, supplied by HealthLink Limited, is accurate and complete. However, as use and interpretation of this document is beyond the control of HealthLink Limited, no liability, either direct or consequential, can be entertained by HealthLink Limited, its agents, or its suppliers.

Document Control

Document Name	Impulse Integration Guide
This Guide is Applicable To	<ul style="list-style-type: none"> - Impulse - HealthLink SIX version 6.6.x
Date/Time Created	13/07/2011 5:00 p.m.
Author(s)	Rajab Nabi

Revision Date	Revision Number	PMS Version	Author(s)	Summary of Changes
15/06/2010	0.1		Rajab Nabi	New guide format
18/11/2010	1.2		Kyle Macdonald	Updated with 6.6.1 RSDAU & LAB2 to PIT option
13/11/2011	1.2		Rajab Nabi	Updated HMS screen Shots.

Table of Contents

Document Control..... 2

Before you begin..... 4

Setting up Healthlink with Impulse 5

Setting up Pathology Laboratory in Impulse..... 8


Before you begin

- 1- Make sure you have installed the latest version of Healthlink SIX Messaging System (HMS version 6.6.x) in your system
- 2- Healthlink needs to be installed on the Impulse server
- 3- If you are unsure about whether you have the latest Healthlink SIX software. Or whether you have Healthlink at all then please contact the Healthlink helpdesk on 180 125 036

Setting up Healthlink with Impulse

Impulse currently accepts reports in PIT format only.


- 1- Open Healthlink advance options by clicking **Start > All Programs > Healthlink SIX > Advance Options**
- 2- Click on **Configuration > User Settings**. The User Settings window will appear.
- 3- Click on **Message Types** tab.
- 4- Scroll down the **Message Specifications** list and select PIT
- 5- Click on the **Directories** tab in the lower part of the screen
- 6- Click on **Edit/More** button. The **Messaging Directories** window will appear.

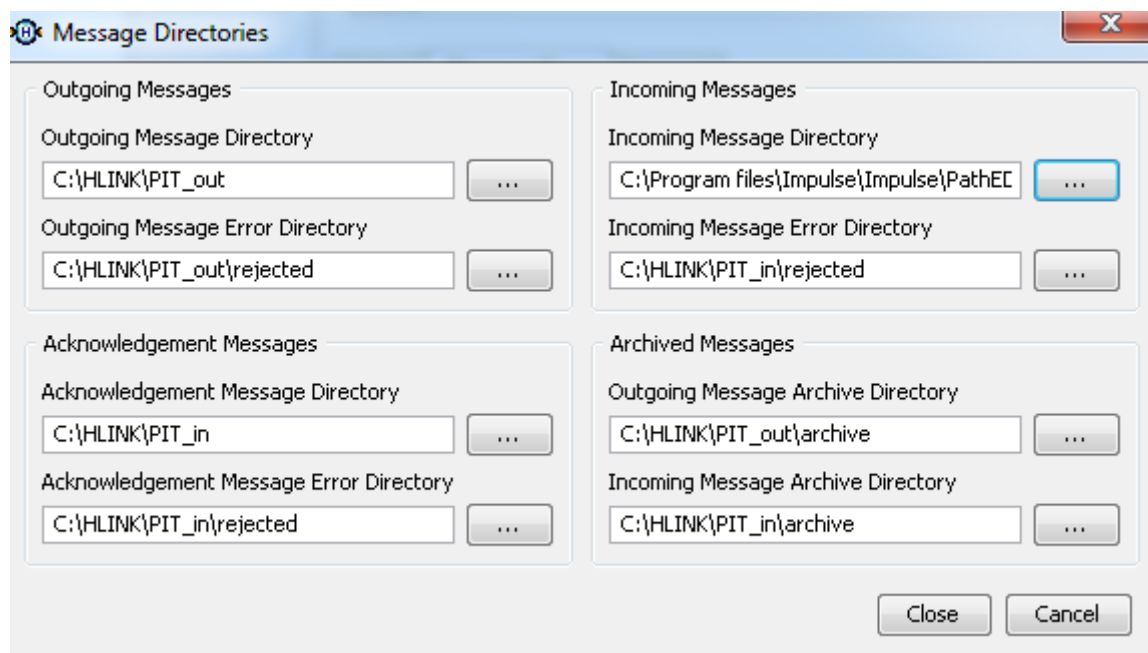

- 7- Change **Incoming Message Directory** to C:\Program files\Impulse\Impulse\PathEDI\Inbound
- 8- Click on **Message Types** tab.
- 9- Scroll down the **Message Specifications** list and select **LAB2**
- 10- Click on the **Directories** tab in the lower part of the screen


11-Click on **Edit/More** button. The **Messaging Directories** window will appear.

12-Update the Outgoing and incoming Message Directories as shown below.

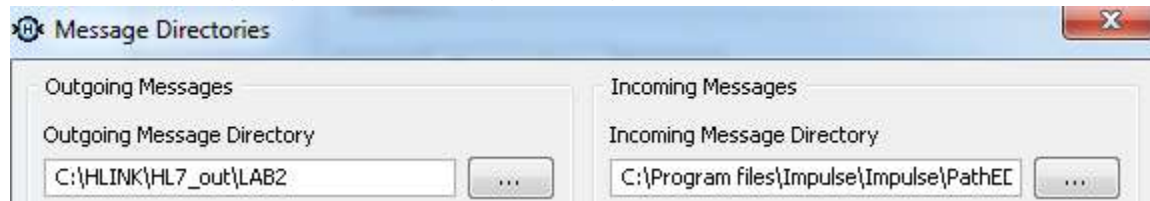

13-Click the **Close** button

14-Click on the **Processor** Tab

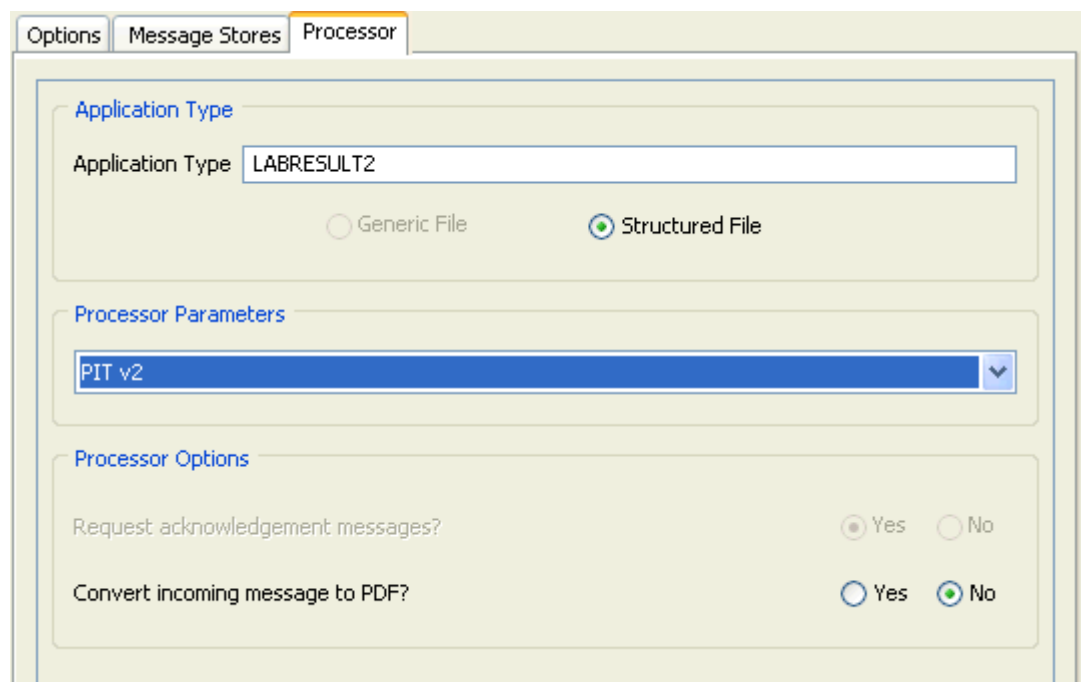

15-In the Processor Parameter choose **PIT v2**.


16-Check that the Convert to PDF is No

17-Scroll down the **Message Specifications** list and select **RSDAU**

18-Click on the **Directories** tab in the lower part of the screen


19-Click on **Edit/More** button. The **Messaging Directories** window will appear.

20-Update the Outgoing and incoming Message Directories as shown below.

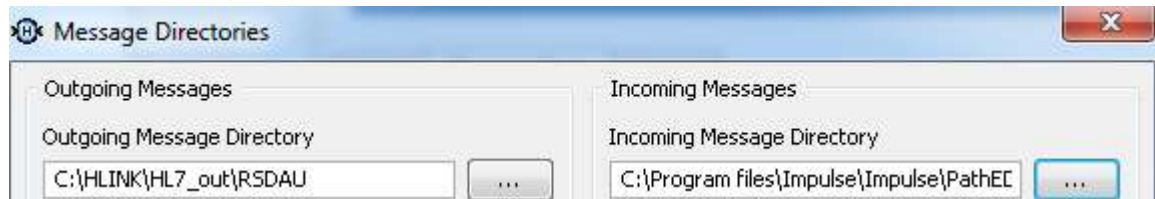

21-Click the **Close** button

22-Click on the **Processor** Tab

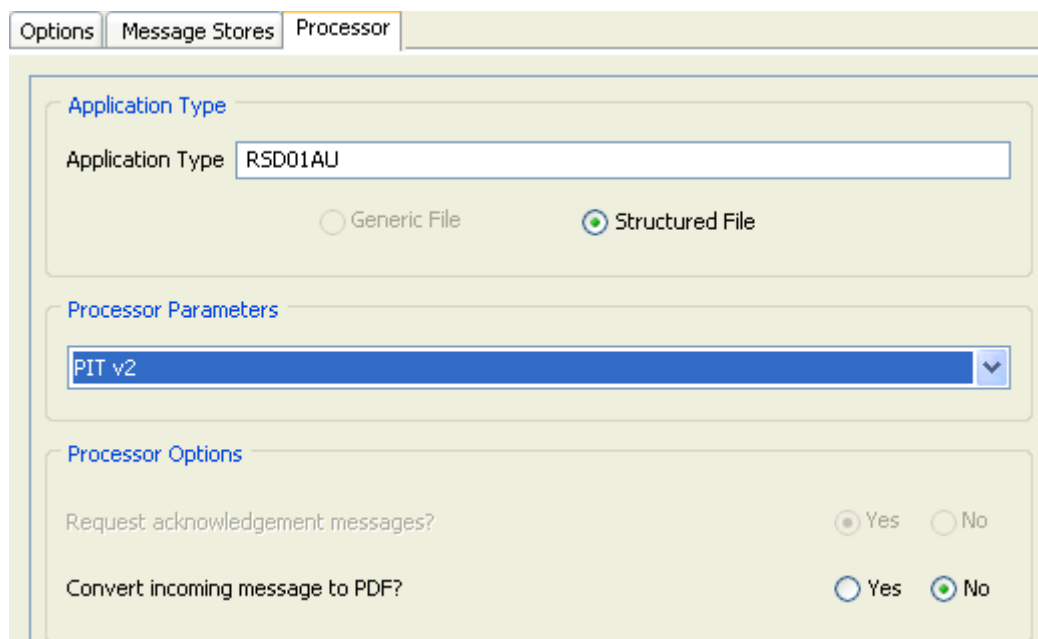

23-In the Processor Parameter choose **PIT v2**.

24-Check that the Convert to PDF is No

25-From the **File** menu, click **Save All**

26-Close Healthlink SIX advance options

Setting up Pathology Laboratory in Impulse

- 1- Click the **Plus** to Add the Lab
- 2- Lab code is usually **3 Alpha 4 Numeric**. Shorten the lab name to initials.
Example *Healthlink Laboratory* will become *HLL00001*
- 3- Copy the first **PIT Line** excluding the line identifier.
Example *Healthlink Pathology Report* *06 15/06/2010*
- 4- Enter the rest of the address details and laboratory type as **Pathology**
- 5- Tick to save

